

LUNA™ Reusable Slide

L12011
L12012

Product Description

Slide

Materials	Aluminum alloy, glass
Dimensions (W x D x H)	75 mm x 25 mm x 2.5 mm
Chamber volume	10 µL


Coverslip

Material	Glass
Dimensions (W x D x H)	20 mm x 15 mm x 0.6 mm

The LUNA™ Reusable Slide is used for cell counting with the LUNA™ family of automated cell counters and CELENA™ S Digital Imaging System. The LUNA™ Reusable Slide may be used for both fluorescence and brightfield cell counting.


Directions for Use

1. Place a clean LUNA™ Reusable Slide Coverslip carefully on top of the LUNA™ Reusable Slide chamber.

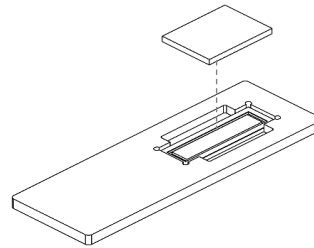

NOTE: Use only the provided coverslip with the LUNA™ Reusable Slide to avoid focusing issues with the LUNA™ family of automated cell counters.

2. Place the pipette into the space between slide and coverslip at the edge of the chamber. Fill with 10-12 µL cell suspension.

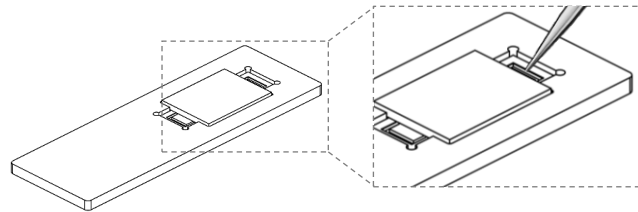

NOTE: Do not lift the coverslip with the pipette tip.

NOTE: Do not overfill.

3. Count the sample with a compatible LUNA™.

Product Care

1. Rinse the LUNA™ Reusable Slide and Coverslip with tap water followed by 70% ethanol. Handle gently to avoid damage.
2. Dry gently and thoroughly with Kimwipes.
3. Wrap in lens paper and store in a safe location.

! Do not press down on the LUNA™ Reusable Slide and Coverslip.
Do not apply excessive force to the slide chamber when handling.

! Do not leave the LUNA™ Reusable Slide and Coverslip in solvents for extended periods of time.

! Do not sonicate the LUNA™ Reusable Slide and Coverslip.

Disclaimer

This product is for research use only.

Additional information is available on our website at www.logosbio.com.
© Logos Biosystems 2018. All rights reserved.



HEADQUARTERS

FL 2 & 3
28 Simindaero 327beon-gil, Dongan-gu
Anyang-si, Gyeonggi-do 14055
South Korea

Tel: +82 31 478 4185

USA

7700 Little River Turnpike STE 207
Annandale, VA 22003
USA

Tel: +1 703 622 4660

EUROPE

11B avenue de l'Harmonie
59650 Villeneuve d'Ascq
France

Tel: +33 (0)3 74 09 44 35

www.logosbio.com

